

Certification



Promoting member of ISMST



CE (MDD)



ISO13485 / ISO9001




Health Canada



Therapeutic Goods Administration



INMETRO

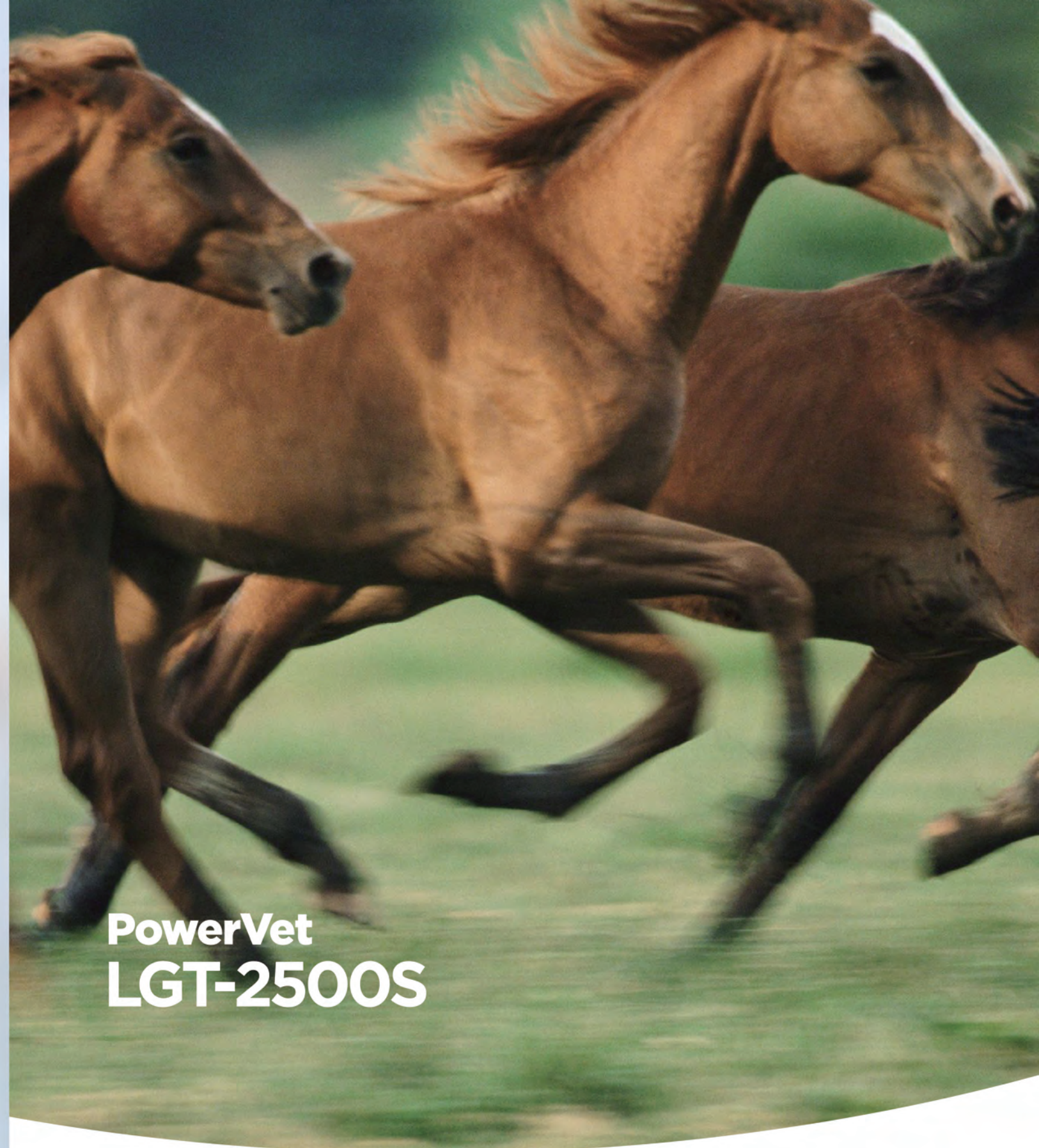

**PowerVet
LGT-2500S**

Longest

Guangzhou Longest Science & Technology Co., Ltd.

Address: 5&6F, Building B4, No.11, Kaiyuan Avenue, Science City, Guangzhou Hi-tech Industrial Development Zone, 510530 Guangzhou, China

Tel: 86-20-66353999 Fax: 86-20-66353920

E-mail: export@longest.cn Website: www.longest.cn

www.longest.cn

Longest

Electrical Extracorporeal Shock Wave Therapy Device (PowerVet LGT-2500S) is the ideal ESWT equipment for equine and small animals use.

The method of shock wave therapy based on short-term application to the field of diseases of high-low-frequency sound wave shock. Extracorporeal shock wave therapy (ESWT) has become an emerging new technology for treating musculoskeletal problems, soft-tissue injuries and bone injuries in horses and small animals. ESWT is non-invasive, used to stimulate healing to return horses and small animals to a level of fully sound, useful activity without recurrence of disease.

Medical studies show that Shock Wave Therapy is effective, non-invasive treatment of orthopaedic and neurological conditions both in race horses and small animals.



PowerVet LGT-2500S




Electrical Extracorporeal Shock Wave Therapy Device



⊕ Technical Specifications

- ⊗ **Pressure:** 1.0 - 4.0 bar (step by 0.1 bar)
- ⊗ **Frequency:** 1 - 17 Hz (step by 1 Hz)
- ⊗ **Compressor type:** Built-in air compressor
- ⊗ **Dimension:** 39.5(L) × 39(W) × 19.4(H) cm
- ⊗ **Weight:** approx 8.5 kg (handpiece excluded)


Radial shock wave handpiece & transmitters

	PowerVet Handpiece ergonomic design durable easy-use
	15R 15mm radial transmitter for all tendinopathies
	20R 20mm radial transmitter for myofascial trigger therapy
	15F 15mm focus transmitter* fascial muscles neck, surface pain area
	15D 15mm deep transmitter* deep-seated disease, trigger point
	6A 6mm acupuncture transmitter* shock wave acupuncture



notice: * stands for "optional accessories"

⊕ Indications for horses:

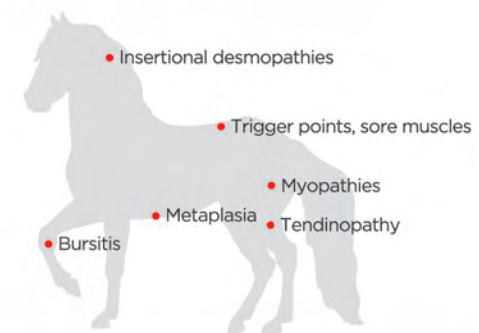

⊕ Indications for dogs:



Make your therapy anywhere
Leave your competition in the dust