

## Certification



National Medical Products Administration



CE (MDD)



ISO13485 / MDSAP

# LGT-2200 Series

## Compression Therapy Device



# Longest

Guangzhou Longest Science & Technology Co., Ltd.

Address: 5&6F, Building B4, No.11, Kaiyuan Avenue, Science City,  
Guangzhou Hi-tech Industrial Development Zone, 510530 Guangzhou, China

Tel: 86-20-66353999 Fax: 86-20-66353920

E-mail: [export@longest.cn](mailto:export@longest.cn) Website: [www.longest.cn](http://www.longest.cn)



Health Life Longest Care


# Longest



**LGT-2200 Series** is sequential pneumatic pressure system that through program controlled inflation and deflation of the sleeves covered in arms, legs, waist, and back to realize therapeutic effect


### ⊕ Mechanism and functions of LGT-2200 Series

The inflation and deflation of the sleeves enhances the backflow of venous blood and lymph, which can improve the blood circulation of diseased areas. This will lead a series of beneficial effect such as improving the function of peripheral blood vessels, removing the edema, accelerating the healing process, and preventing Muscular Dystrophy.

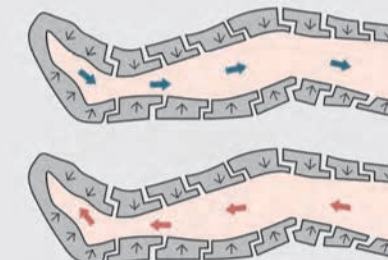

### ⊕ Features

- ⊗ Overlapping air chambers generate a gliding pressure wave with high efficiency
- ⊗ User-friendly interface, just a few steps to finish your setting
- ⊗ Gentle inflation and deflation
- ⊗ Real-time pressure, time, and inflation part monitoring (visible in the interface)
- ⊗ Durable & easy-use sleeves
- ⊗ Pressure unit switchable between mmHg and kPa

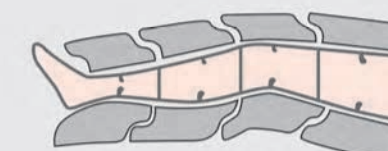

Overlapping Air Chamber



### ⊕ Specifications



Professional

#### LGT-2200L

- ⊗ **Chambers:** 8
- ⊗ **Compression modes:** 18 ( Including 3 DVT modes + 6 four chambers modes + 9 eight chambers modes )
- ⊗ **Time range:** 1 - 99 mins
- ⊗ **Pressure range:** 40 - 190 mmHg
- ⊗ **Interface:** 5.7" touch screen
- ⊗ **Dimension:** 393 (L) x 366 (W) x 194 (H) mm
- ⊗ **Weight:** 10 kg
- ⊗ **Input voltage:** AC 220 - 240 V, 50 Hz ( LGT-2200L)  
AC 100 -120 V, 60 Hz ( LGT-2200LL)



Portable

#### LGT-2200S

- ⊗ **Chambers:** 4
- ⊗ **Compression modes:** 6
- ⊗ **Time range:** 1 - 99 mins
- ⊗ **Pressure range:** 40 - 190 mmHg
- ⊗ **Interface:** 5.7" touch screen
- ⊗ **Dimension:** 315 (L) x 215 (W) x 251 (H) mm
- ⊗ **Weight:** 6.0 kg
- ⊗ **Input voltage:** AC 220 - 240 V, 50 Hz ( LGT-2200S)  
AC 100 - 120 V, 60 Hz ( LGT-2200SL)

### ⊕ Application

- 1/Prophylaxis of Deep vein thrombosis(DVT) (Post-stroke deep vein thrombosis 、 Post-traumatic deep vein thrombosis)
- 2/Prophylaxis and treatment of edema (Chronic venous edema、 Lymphedema 、 Post-mastectomy lymphedema)

