

## POWER-WAVE LGT-2510 Series

Radial Acoustic Wave Therapy


### Acoustic Wave Therapy

Acoustic Wave Therapy is a type of treatment that applies mechanical vibrations to the outside surface of the skin. These pulsations stimulate the treated area in a non-invasive manner, meaning they can penetrate through the surface of the skin without breaking it. As an extension application of Shock Wave Therapy in Aesthetic medicine, Acoustic Wave Therapy has been widely used in the treatment of cellulite and other skin irregularity problems such as stretch marks and unsightly scars. After years of clinical experiment, it has been proved that Acoustic Wave Therapy is a safe and effective therapy. LONGEST POWER-WAVE LGT-2510 series takes advanced pneumatic ballistic technology and use high-performance oil compressor, which guarantees stable powersource supply and superior treatment effect. It has also a wide range of application.

### Application of POWER-WAVE LGT-2510 Series

- Anti-cellulite treatment
- Body reshaping
- Skin elasticity improvement
- Connective tissue tightening
- Scar and wrinkle smoothing
- Treatment of stretch marks
- Lymphatic system stimulation
- Treatment of skin irregularity after liposuction

# EASY TO REBUILD YOUR BEAUTIFUL LIFE! POWER-WAVE LGT-2510 SERIES

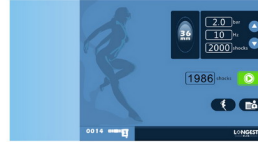


#### HIGH ACCURATE PARAMETER ADJUSTMENT

- Parameters adjustable during the treatment
- Pressure steps low to 0.1 bar
- Frequency steps low to 1.0 Hz




#### PRE-SET PROTOCOLS

- Imaged-based treatment guidance
- Pre-set protocols based on clinical tests
- One-touch-start treatment

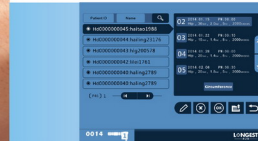

#### PATIENT DATABASE WITH CIRCUMFERENCE

- Up to 5000 patient files to be stored
- One touch search of the treatment record
- Circumference record for easy tracking of treatment effect

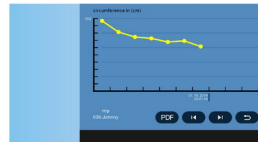

#### ALL SESSION AT A GLANCE

- Circumference record in patient file
- Curve diagram to compare treatment effect
- PDF file output supported

Model Number:  
LGT-2510A

#### Advantages

- Non-surgical
- Non-invasive
- No side-effects
- Anesthesia-free
- Pain-free



Upper Arm



Thigh




Abdomen



Hip



Buttocks

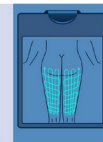

Thigh

## POWER-WAVE LGT-2510 Series




Radial Acoustic Wave Therapy

### LONGEST Designed Handpiece

- New design
- Ergonomic
- Durable & easy-use
- Soft & comfortable feeling



### Transmitter

	Diameter	Function
	D=36mm Radial	Smoothing muscles and connective tissue over large areas; stimulating blood and lymph flow.
	D=20mm Radial	Stimulating blood and lymph flow; enhancing collagen.
	D=15mm Deep	Breaking down fat and connected tissues in deep-seated areas.

## TECHNICAL SPECIFICATIONS

Model Number	LGT-2510A	LGT-2510B
Channels	1	2 independent
Interface	10" LCD touch screen	10" LCD touch screen
Pressure range	1.0~5.0 bar	1.0~5.0 bar
Frequency range	1~22Hz	1~22Hz
Pre-set protocols	Yes	Yes
Patient database	Yes	Yes
Continuous shock mode	Yes	Yes
Special pulsed mode	Yes	Yes
Software update	Yes	Yes
Item dimension	54.2*44.5*112cm	54.2*44.5*112cm
Item weight	58kg	58kg
Input voltage	110~120V / 60 Hz, 220~240V / 50 Hz	
Safety standard	Class II a	Class II a

## ACCESSORIES

Model Number	LGT-2510A	LGT-2510B
Handpiece	1	2
36mm Radial Transmitter	1	1
20mm Radial Transmitter	1	2
15mm Deep Transmitter	1	2
Projectile	2	4
Metal Tube	2	4
Gel 250 ml	1	1
Reversion Box	1	2



Promoting Member of ISMST

## Power-Wave LGT-2510 Series

Radial Acoustic Wave Therapy



GUANGZHOU LONGEST SCIENCE&TECHNOLOGY CO., LTD

Address: 5&6F, Building B4, No. 11, Kaiyuan Avenue, Science City  
Guangzhou Hi-tech Industrial Development Zone, 510530 Guangzhou  
Office Tel: +86) 020 6635 3999 Email: export@longest.cn  
Office Fax: +86) 020 6635 3920 Website: www.longest.cn



Authorized distributor: